




## Measuring your dog for a set of Walkin' Wheels

We only need your dog's weight, breed and one height measurement.

1. With your dog supported in their correct standing position, spine parallel to the ground, measure from the ground to the groin. This is the fold in the flank of the dog, the web like skin that joins the back leg to the abdomen (look at the skin between your thumb and forefinger).

If your pet is unable to stand, measure as shown below. Do not straighten your pet's leg. Leave the natural curve.

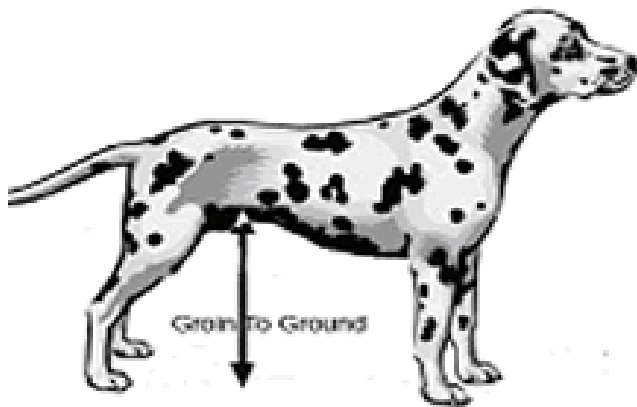

**\*Note:** If your pet is overweight, then the fold of flank will probably hang down, making the fold of flank measurement shorter than it really is. You will need to try and move the flank up to a normal position, to get an accurate measurement. The picture of the Boxer shows a normal fold of flank for a dog that is not overweight.

If you are unsure about how to measure, please call us (0800 072 3697) and we will be happy to help you.

**All you need to do is take the measurement and ask for a quote**